

Hipper — Therapist Documentation

Team Digital Life - Hipper Monitoring

Michel Oey
Indigo van Hoof
Jan Vos
Mike Schardijn
Bryan van Riesen
Richard van Breugel
Victor Hoegy
Darian Kardjo

DIGITAL
LIFE



Table of Contents

1. Introduction	4
1.1 Purpose of the therapist Interface	4
1.2 Who Should Use This Guide	4
2. Getting Started	5
2.1 therapist Access and Login	5
2.2 Navigation Overview	6
3. Managing Patients	7
3.1 Viewing Assigned Patients	7
3.2 Adding a New Patient	7
3.3 Resetting Patient Passwords	8
4. Managing Devices	9
4.1 What is the Movement Score?	9
4.2 Adding or Removing Devices	10
4.3 Troubleshooting Device Issues	10
5. Therapist Preferences	11
5.1 Interface Personalization	11
5.2 Saving Preferences	11
6. Data and Security	14
6.1 What Data Therapists Can See	14
6.2 Privacy Considerations	14
7. Troubleshooting	15
7.1 Common therapist-Side Issues	15
7.2 When to Contact Support	15
8. Contact and Support	17
8.1 Helpdesk Information	17
8.2 Support Hours, Expectations and feedback	17

1. Introduction

1.1 Purpose of the therapist Interface

The Hipper therapist Interface is designed to help therapists and therapististrative staff monitor and support patient recovery through activity tracking. This guide explains how to use the therapist tools to manage users, link devices, and adjust preferences.

1.2 Who Should Use This Guide

The primary target audience for Hipper includes individuals who are recovering from physical This guide is intended for:

- Therapists working with patients using Hipper
- therapististrative users who help manage accounts
- Superusers with elevated system privileges


No technical knowledge is required to use the therapist dashboard.

2. Getting Started

2.1 therapist Access and Login

To access the therapist panel:

1. Go to the Hipper therapist login page.
2. Enter your email and password.
3. After successful login, you will be directed to the therapist dashboard.

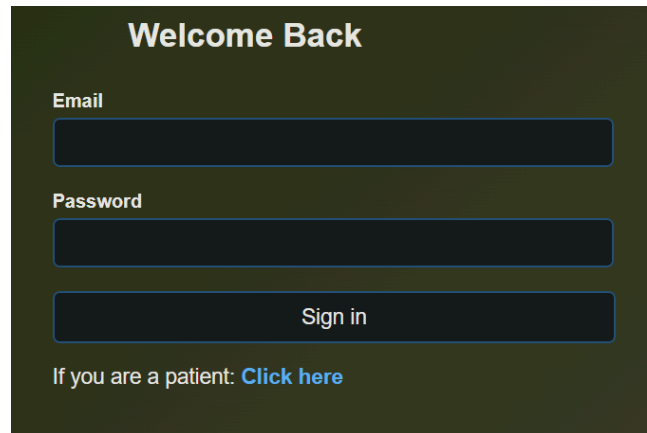

Welcome Back

Email

Password

Sign in

If you are a patient: [Click here](#)

If you forget your password, contact your system therapististrator or platform manager.

3. Managing Patients

3.1 Viewing Assigned Patients

All patients assigned to you are listed in the **Patients** section.

Use the search bar to quickly locate individuals by name.

Each patient entry shows:

- Name
- Email address




The screenshot shows a dark-themed interface titled "Patients". At the top, there is a search bar with the placeholder text "Search patients by name...". Below the search bar, a patient entry is displayed in a rounded rectangle. The entry shows the name "Henk Man" in bold, followed by the email address "Email: henk.man@gmail.com".

3.2 Adding a New Patient

To add a new patient:

1. Go to the **Add Patient** form.
2. Enter the patient's full name, email address, and a secure password.
3. Click **Add Patient**.

If the email is already in use or there is a form error, a red message will appear below the button.

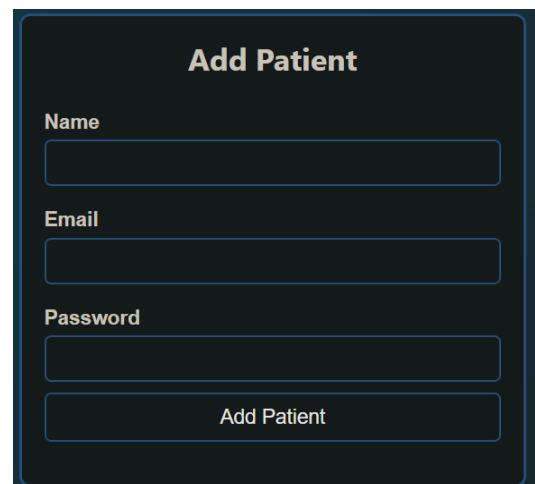


The screenshot shows a dark-themed interface titled "Add Patient". It contains three input fields: "Name", "Email", and "Password". Below the "Password" field is a button labeled "Add Patient".

3.3 Resetting Patient Passwords

If a patient forgets their password:

1. Go to the **Change Patient Password** section.
2. Enter their email and a new password.
3. Click **Change Password**.

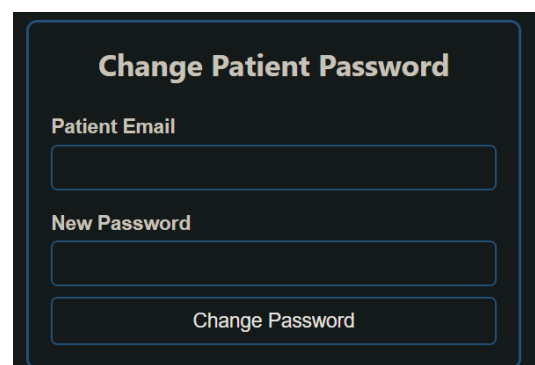

The screenshot shows a dark-themed interface titled "Change Patient Password". It contains two input fields: "Patient Email" and "New Password". Below the "New Password" field is a button labeled "Change Password".

This update takes effect immediately.

4. Managing Devices

4.1 Viewing and Assigning Devices

In the **Manage Devices** section, you can:

- View all active PAM devices
- See which patient each device is assigned to
- Filter by “unlinked” devices for setup

4.2 Adding or Removing Devices

A higher score typically means more movement or more intense movement. A lower score suggests minimal activity. It's important to interpret the score over time rather than based on a single moment—gradual improvement is key.

To add a new PAM device:

1. Enter the label and MAC address.
2. Click **Add Device**.

To assign a device to a patient:

1. Choose a device and a patient from the dropdowns.
2. Click **Assign**.

To unlink or delete a device:

- Use the buttons next to each device entry.

Deleting a device is permanent. Only do this if it's no longer used.

The screenshot displays a web interface for managing PAM devices, organized into three main panels:

- Patients Panel:** Features a search bar labeled "Search patients by name...". Below it, a card for "Patient ID: 1" displays the following information: Name: Henk Man, Email: henk.man@gmail.com, and Bound to Device: 09234.
- (Un)bind device to patient Panel:** Contains two input fields for "Patient ID" and "Device ID". Below these fields are two buttons: "Bind" and "Unbind".
- Devices Panel:** Includes a search bar labeled "Search devices by label...". A button labeled "Show Only Unoccupied Devices" is present. Below this, a card for "Device ID: 1" shows a "Delete" link, "Device Label: 09234", and "Device address: C1:08:00:01:23:80".

4.3 Troubleshooting Device Issues

If a device is not syncing:

- Make sure it is powered on and worn correctly.
- Check the base station (must be within Bluetooth range).
- Use the unassign/reassign process to force a reset.
- If issues persist, contact technical support.

5. Therapist Preferences

5.1 Interface Personalization

Each therapist can adjust the interface based on personal comfort.


Options available:

- **Theme:** Light / Dark
- **Font Size:** Normal / Big

5.2 Saving Preferences

After selecting your preferences, click **Save Changes**.

Preferences are saved across sessions and devices using your account settings.

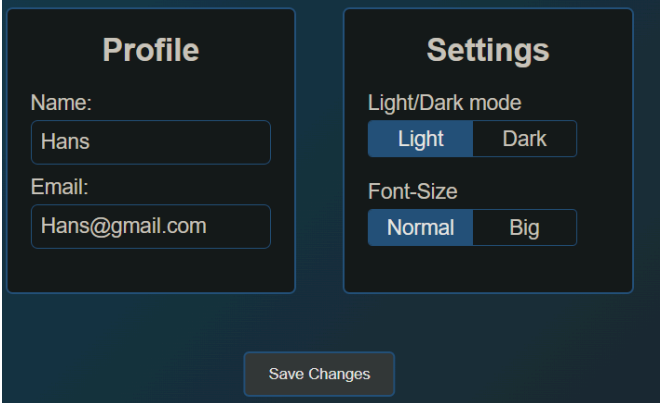


The screenshot displays a user interface for therapist preferences. It is divided into two main sections: 'Profile' and 'Settings'. The 'Profile' section contains two text input fields: 'Name:' with the value 'Hans' and 'Email:' with the value 'Hans@gmail.com'. The 'Settings' section contains two toggle-like options: 'Light/Dark mode' with 'Light' selected and 'Dark' as an alternative, and 'Font-Size' with 'Normal' selected and 'Big' as an alternative. At the bottom center of the interface is a button labeled 'Save Changes'.

6. Data and Security

Therapists only have access to patients assigned to them.

They can view:

- Movement history
- PAM scores and steps
- Device sync status
- Goal progression

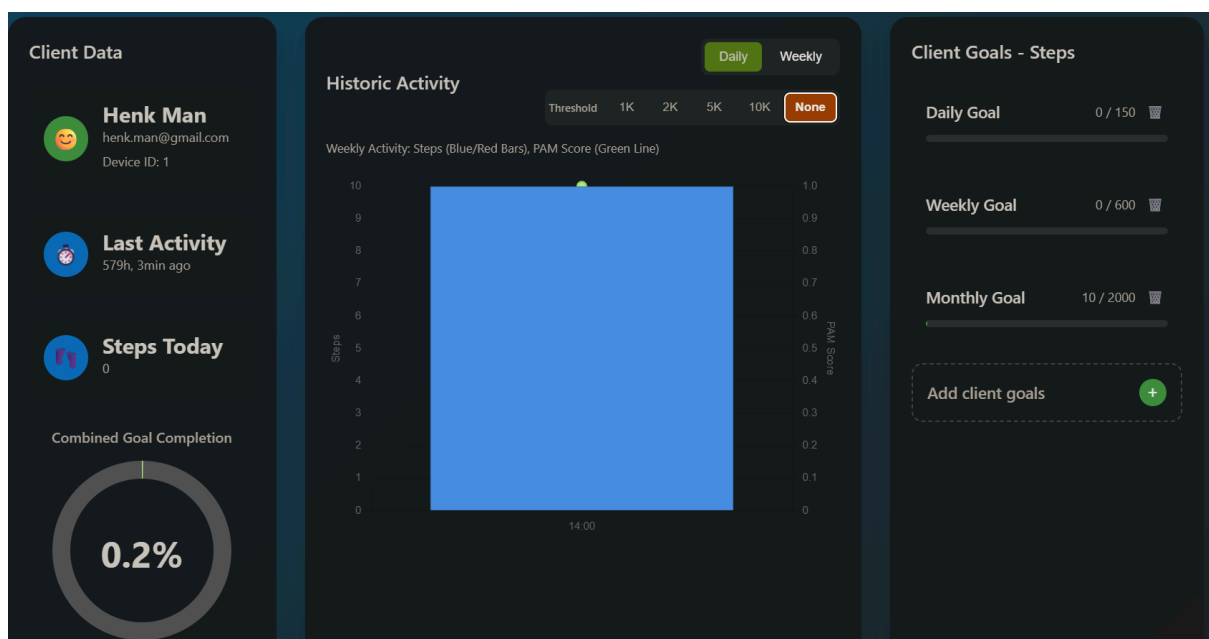

6.2 Privacy Considerations

Hipper complies with data protection laws. Patient data is:

- Only visible to assigned therapists
- Stored securely
- Not shared with third parties

Superusers have limited therapist administrative oversight for compliance purposes.

- This ensures that your activity data remains private, especially when using shared or public devices.

6.3 Super Users

Super users are a special type of user that can manage therapists. instead of seeing “their own patients” they will instead be able to see ALL the patients and ALL the pam devices and manage these patients and pam devices.

but most importantly, super users can add new therapist users and manage them with actions like changing credentials and promoting. since the super users are the highest form of authority among the accounts, they can also add or remove other super users.

Instead of the normal profile settings (which are still available at the bottom of the screen) the super user gets these new functionalities in order to manage therapist and super users.

————— HIPPER THERAPEUTICS —————

You are a super user, This means that you can manage all other super users and therapists!

Manage Super Users:

Promote therapist to super user:

Enter therapist's email...

Active Super Users:

super

super@gmail.com

Demote

test

test@gmail.com

Demote

Manage Therapists:

Add New Therapist

Name: Email:

Password:

Search for therapist by name or email:

Search therapists...

All Therapists:

hans

hans@gmail.com

super

super@gmail.com

test

test@gmail.com

Reset Therapist Password

Therapist's email...

New password...

Profile

Name:

Email:

Settings

Light/Dark mode:

Font-Size:

6.3.1 Manage Super Users

Promote Therapist to Super User

- Enter the therapist's email address in the "Promote therapist to super user" field.
- Click the "Promote" button.
- The therapist will be added to the Active Super Users list.

Active Super Users

- The list shows all current Super Users and their email addresses.
- To revoke Super User status, click "Demote" next to the user.

6.3.2 Manage Therapists

Add New Therapist

- Complete the Name, Email, and Password fields under "Add New Therapist."
- Click "Add Therapist" to register the new therapist.

Search for a Therapist

- Use the search box under "Search for therapist by name or email."
- Typing a name or email filters the list in real time.

All Therapists

- Displays every registered therapist and their email.
- To remove a therapist, click "Delete" next to their entry.

Reset Therapist Password

- Enter the therapist's email in the "Therapist's email..." field.
- Enter a new password in the "New password..." field.
- Click "Reset Password" to update their credentials.

7. Troubleshooting

7.1 Common Issues

1.

Problem	Solution
Patient not appearing	Check if they are assigned to you
Device won't assign	Ensure MAC is correct and not already linked
Preferences won't save	Refresh the page and try again

Dashboard doesn't load

Clear cache and check internet connection

7.2 When to Contact Support

Contact support if:

- You can't log in or reset your password
- Devices stop syncing completely
- Patients are missing or corrupted in the system

8. Contact and Support

8.1 Helpdesk Information

If you encounter any issues with the dashboard, device syncing, or your account, our helpdesk is here to assist you.

- Email: support@hipperacademy.nl

8.2 Support Hours

Our team is available on weekdays from 09:00 to 17:00 (CET).

Outside these hours, you're welcome to leave a message — we'll get back to you as soon as possible.

8.3 Response Time

We aim to respond to most support requests within 1–2 business days.

Critical issues, such as login problems or missing data, are prioritized and typically addressed within 24 hours.

9.4 Feedback Submission

Your feedback is always appreciated. If you'd like to report a bug, make a suggestion, or help us improve the platform:

- Send us a message at: info@hippertext.nl