

Hipper — User Documentation

Team Digital Life - Hipper Monitoring

Michel Oey
Indigo van Hoof
Jan Vos
Mike Schardijn
Bryan van Riesen
Richard van Breugel
Victor Hoegy
Darian Kardjo

DIGITAL
LIFE



Table of Contents

1. Introduction	4
1.1 Purpose of the project	4
1.2 Target audience	4
1.3	
2. Getting Started	5
2.1 System Requirements	5
2.1.1 Supported Browsers	5
2.1.2 Internet Connection	5
2.1.3 Screen Resolution Recommendations	5
2.2 Creating an Account	5
2.3 Logging In	5
2.4 Accessing the Dashboard	6
3. Dashboard Overview	7
3.1 Complete Overview	7
3.2 Today's Activity	7
3.3 Combined Goal Completion	8
3.4 Personal Goals — Steps	8
3.5 Historic Activity	8
4. Movement Score Explained	9
4.1 What is the Movement Score?	9
4.1.1 How it's Calculated	9
4.1.2 Units of Measurement	9
4.1.3 Normal Score Ranges	9
4.2 How to Interpret the Score	9
4.2.1 What High or Low Scores Mean	9
4.2.2 Color Coding	10
4.2.3 Daily Goals or Thresholds	10
4.3 Progress Tracking	10
4.3.1 Reading Charts and Trends	10
4.3.2 Filtering by Date or Time	10
5. Device Integration	11
5.1 How Data is Collected	11
5.1.1 Description of the Waist Device	11
5.2 Synchronizing of data	11
5.2.1 How and When It Syncs	11
5.2.2 Troubleshooting Syncing Issues	11
6. Account Management	12
6.1 Profile Settings	12
6.1.1 Updating personal details	12
6.1.2 Managing preferences	12
6.2 Security	12

6.2.1 Changing passwords	12
6.2.2 Logging out	12
6.3 Deleting account	12
6.3.1 How to request account deletion	12
6.3.2 Data retention and privacy policies	13
7. Privacy and Data Security	14
7.1 What Data is Collected	14
7.2 How Data is Used	14
7.3 Your Rights	14
7.3.1 Data Access, Correction, and Deletion	14
8. Troubleshooting	15
8.1 Common Issues	15
8.1.1 Dashboard Not Loading	15
8.1.2 Missing Data	15
8.1.3 Inaccurate Scores	15
8.2 Device-Related Issues	15
8.2.1 Device Not Syncing	15
8.2.2 Battery Issues	16
8.3 Who to Contact If the Issue Persists	16
9. Contact and Support	17
9.1 Helpdesk Contact	17
9.2 Support Hours	17
9.3 Response Time	17
9.4 Feedback Submission	17

1. Introduction

1.1 Purpose of the project

Hipper is a project developed and maintained by Digital Life, aimed at creating a movement tracker for individuals undergoing rehabilitation or for the elderly. The goal of this project is to offer a straightforward yet reliable way for users to monitor their daily physical activity. By making it easier to stay aware of their movements throughout the day, Hipper hopes to support users in maintaining an active lifestyle, which in turn can contribute to faster recovery and improved overall well-being.

1.2 Target audience

The primary target audience for Hipper includes individuals who are recovering from physical injuries, patients undergoing rehabilitation, and elderly people who benefit from increased awareness of their daily movement. These groups often face challenges in maintaining consistent levels of activity, either due to limited mobility, health concerns, or simply the need for structured support.

Hipper is designed with these users in mind—offering a tool that is both accessible and easy to use, without requiring technical know-how. Whether it's helping someone rebuild strength after surgery or encouraging gentle movement to maintain independence in older age, Hipper is tailored to meet the needs of people who can gain the most from careful and consistent activity tracking.

2. Getting Started

2.1 System Requirements

2.1.1 Supported Browsers

Hipper is compatible with all modern web browsers, including Google Chrome, Safari, Microsoft Edge, and the latest versions of Firefox. For the best experience, we recommend keeping your browser up to date to ensure full functionality and security.

2.1.2 Internet Connection

Hipper requires an active internet connection to function properly. A stable Wi-Fi or mobile data connection is recommended to ensure real-time syncing of movement data and access to your personal activity dashboard. The platform is designed to handle typical household internet speeds without issue.

2.1.3 Screen Resolution Recommendations


While Hipper can be accessed on various screen sizes, we recommend using a screen with a minimum resolution of 1280 x 720 pixels for optimal readability and layout. The interface is responsive and works well on tablets, laptops, and desktop displays. Mobile phones are supported, but for the best user experience, a larger screen is ideal when viewing detailed activity summaries.

2.2 Creating an Account

Patients are not required to create an account themselves. Instead, their therapist will set up the account on their behalf using the patient's email address and a secure password. This ensures a smooth onboarding process and helps maintain consistency across patient records.

2.3 Logging In

Once your therapist has created an account for you, logging in is simple. Go to the Hipper login page and enter the email address and password provided by your therapist. After logging in, you'll be taken directly to your personal dashboard, where you can view your movement data and track your progress.

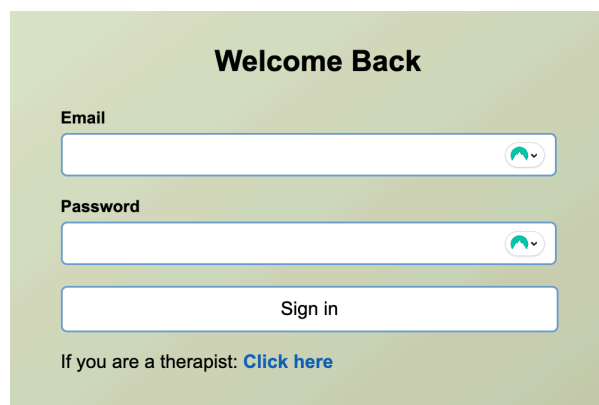

The screenshot shows a login interface with a light green background. At the top, the text "Welcome Back" is displayed in bold. Below this, there are two input fields: "Email" and "Password". Each field has a small green icon with a white arrow pointing right. Below the password field is a "Sign in" button. At the bottom, there is a link that says "If you are a therapist: [Click here](#)".

If you forget your password, please contact your therapist for a password reset.

2.4 Accessing the Dashboard


After successfully logging in, you'll be taken straight to your personal dashboard. This is the central place where you can view all your movement data in a clear and organized way. The dashboard shows your daily activity levels, historical trends, and any relevant notes or feedback from your therapist.

Navigation is designed to be simple and user-friendly, so you can focus on tracking your progress without distractions. From here, you can also access settings, view past activity logs, or log out when you're finished.

3. Dashboard Overview


3.1 Complete Overview

The dashboard is divided into four main areas, each designed to give the patient a clear and meaningful overview of their physical activity. Together, these sections help patients stay informed, motivated, and on track with their recovery or fitness goals.

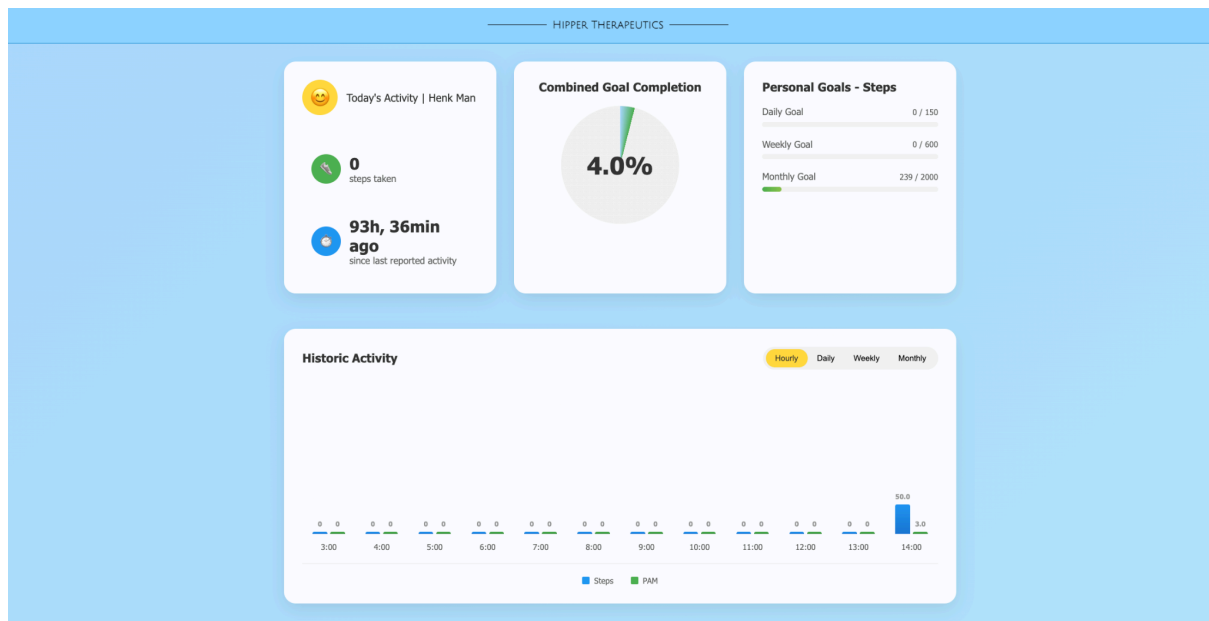


3.2 Today's Activity

This section shows a real-time summary of your movement for the current day. It includes step count and the inactivity in minutes.

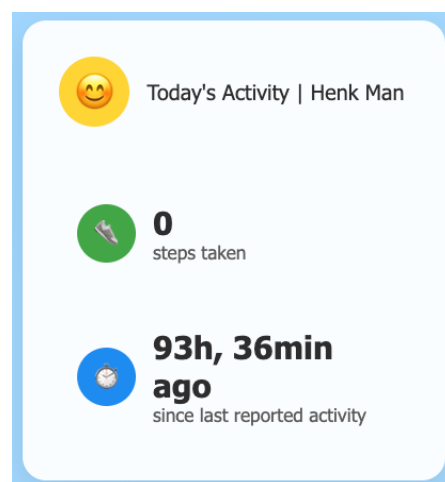

3.3 Combined Goal Completion

Here you'll find a visual representation of how close you are to meeting your daily activity goals. Whether its steps taken or minutes moved, this section helps keep you motivated by showing your progress as it happens.

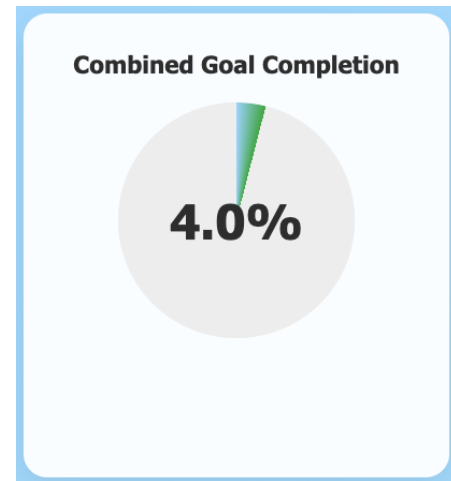


3.4 Personal Goals — Steps

In this area, you can view the activity goals set specifically for you—either by your therapist or in collaboration with them. These goals are tailored to your personal situation and rehabilitation plan, making them realistic and meaningful.


The goals present are:

- Daily goal
- Weekly goal
- Monthly goal

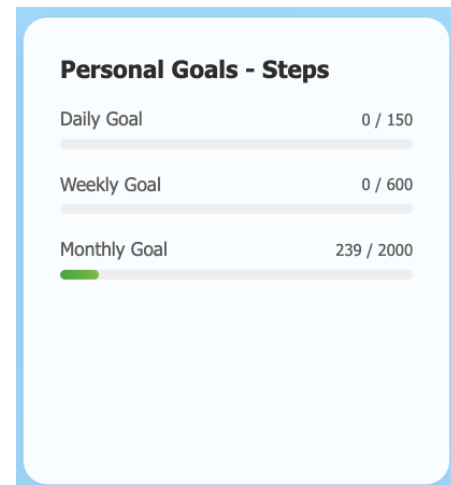

3.5 Historic Activity

This section lets you look back at your movement data over time. You can review past days, weeks, or even months, helping you and your therapist spot trends, improvements, or areas where additional support might be needed. The blue bar is the steps and the green bar is the PAM score accompanied by those steps.

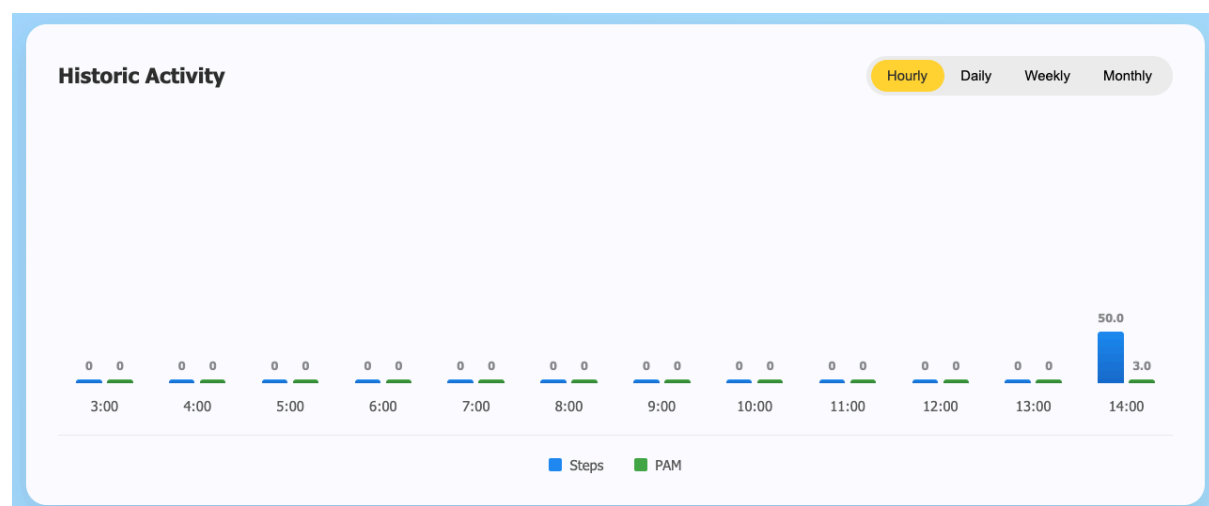

4. Movement Score Explained

4.1 What is the Movement Score?

The Movement Score—often shown as PAM on the dashboard—is a numerical value that reflects how physically active a person has been over a given period. It gives more detail than just step count by taking into account movement intensity and duration.

4.1.1 How it's Calculated

The Movement Score is calculated using data from a wearable device, combining the amount and intensity of movement. It may factor in how long you've been moving, how vigorous that movement was, and how consistently it occurs throughout the day.

4.1.2 Units of Measurement

The score is displayed as a unit-less number (e.g., PAM: 3.0) which increases with more activity. Unlike steps, it's not tied to a specific physical measurement, but is a calculated value based on a custom algorithm.

4.1.3 Normal Score Ranges

While specific goals may vary per patient, generally:

- 0–10 PAM indicates low activity
- 10–50 PAM indicates moderate activity
- 50+ PAM suggests a high level of movement

Your therapist may provide personalized targets within or outside this range depending on your condition.

4.2 How to Interpret the Score

A higher score typically means more movement or more intense movement. A lower score suggests minimal activity. It's important to interpret the score over time rather than based on a single moment—gradual improvement is key.

4.2.1 What High or Low Scores Mean

- High scores suggest you're staying active, which supports recovery or maintenance.
- Low scores may mean you've been inactive and could benefit from gentle movement or consultation with your therapist.

4.2.2 Color Coding

In the dashboard:

- Green bars represent PAM scores (movement intensity)
- Blue bars represent step count

Color may also change based on goal completion (e.g., green if you're meeting your goal).

4.2.3 Daily Goals or Thresholds

Your daily PAM goal is set by your therapist. This will vary depending on your rehabilitation plan. Meeting your daily goal is shown visually in the "Combined Goal Completion" chart and the "Personal Goals" section.

4.3 Progress Tracking

The dashboard provides real-time updates on your Movement Score and allows you to see how you're progressing throughout the day. This keeps you informed and motivated.

4.3.1 Reading Charts and Trends

In the Historic Activity section, you can track your movement:

- Hourly: See what time of day you're most active
- Daily/Weekly/Monthly: Monitor long-term trends

Each bar shows both Steps (blue) and PAM (green) for easy comparison.

4.3.2 Filtering by Date or Time

Use the filter buttons at the top right of the Historic Activity chart to switch between Hourly, Daily, Weekly, or Monthly views. This helps you zoom in on specific timeframes and understand patterns over time.

5. Device Integration

5.1 How Data is Collected

Movement data is collected using a wearable waist device that continuously tracks physical activity throughout the day. This device is designed to be lightweight, discreet, and comfortable enough to wear during regular daily routines.

5.1.1 Description of the Waist Device

The waist device monitors movement patterns and intensity, translating them into step counts and Movement Scores (PAM). It is worn securely on the waist to ensure accurate tracking of core body movement, which is a reliable indicator of overall physical activity.

5.2 Synchronizing of data

5.2.1 How and When It Syncs

The device connects via Bluetooth to a nearby base station, typically located in the patient's home or care facility.

- The base station collects data from the waist device once per hour.
- After collecting the data, the base station uploads it automatically to the server.
- From there, the data becomes visible on the patient's dashboard.

This process is automatic and requires no interaction from the patient.

5.2.2 Troubleshooting Syncing Issues

If you notice that your data hasn't updated in several hours or days, try the following:

- Make sure the device is charged and worn properly on the waist.
- Confirm that the base station is powered on and within range (ideally within 5–10 meters).
- Avoid placing the base station near thick walls or metal objects that may interfere with Bluetooth signal.
- If possible, move closer to the base station and wait 10–15 minutes to allow syncing.

If the problem persists, contact your therapist or support team for further assistance.

6. Account Management

6.1 Profile Settings

Each patient account includes basic profile settings that may be configured by the therapist at setup. These typically include your name, and email address.

6.1.1 Updating personal details

If your personal details change (e.g., email address or name), updates should be made through your therapist. Patients do not have direct access to change these details themselves to ensure data consistency and security. Your therapist can make adjustments quickly through their admin interface.

6.1.2 Managing preferences

Any custom preferences—such as theme, font size—can be changed under the settings tab. These help tailor your experience to match your specific needs and goals.

6.2 Security

Your account is protected by a secure login system using your email address and a password.

6.2.1 Changing passwords

Changing passwords can be done by contacting your therapist, they have the option to change your password.

6.2.2 Logging out

To log out:

- Simply click the “Log Out” button, usually located in the upper-right corner of the dashboard.
- This ensures that your activity data remains private, especially when using shared or public devices.

6.3 Deleting account

Patients do not have the ability to delete their own account directly, to avoid accidental loss of important rehabilitation data.

6.3.1 How to request account deletion

If you wish to delete your account, contact your therapist or support staff. They will guide you through the process and verify your request for security purposes.

6.3.2 Data retention and privacy policies

Hipper follows strict privacy and data protection regulations.

- Your data is only used for medical and therapeutic purposes, and is never shared with third parties without consent.
- Even after account deletion, some data may be retained temporarily as required by law or medical record-keeping standards.
- You can request a full copy of your data at any time, or ask for a detailed explanation of what is stored and why.

7. Privacy and Data Security

7.1 What Data is Collected

Hipper collects a limited set of data necessary to support personal recovery and activity tracking. This includes:

- Personal Information: Your name and email address, provided by your therapist.
- Movement Data: Step counts, Movement Score (PAM), timestamps, and device sync logs.
- System Metadata: Such as the last sync time, device status, and platform usage patterns.

No sensitive health data (such as diagnoses or medication) is stored in the system unless explicitly added by a therapist.

Types of Personal and Movement Data

- Personal Data: Email address, display name, and account ID.
- Activity Data: Daily steps, PAM score, goal progress, and historic activity logs.
- Session Data: Device synchronization time.

7.2 How Data is Used

All collected data serves the following purposes:

- Personal Monitoring: To help you and your therapist track activity progress over time.
- Rehabilitation Support: Allows therapists to set goals and adjust recovery plans based on accurate movement data.

7.3 Your Rights

As a user, you have full control over your personal data. This includes the right to:

- Access your data and view how it's used.
- Request corrections to any inaccurate information.
- Request deletion of your account and associated data (see Account Management for how).
- Withdraw consent for optional data uses, such as anonymized research.

7.3.1 Data Access, Correction, and Deletion

If you'd like to review, change, or delete your data:

- Contact your therapist or platform administrator.
- Requests are typically processed within 14 days.
- In some cases, limited data retention may apply for legal or medical reasons.

8. Troubleshooting

This section outlines common issues users may encounter while using Hipper and provides practical steps to resolve them. If problems persist after following the suggested solutions, you should reach out to your therapist or support contact.

8.1 Common Issues

8.1.1 Dashboard Not Loading

Cause: Poor internet connection, outdated browser, or temporary server issue.

Solution:

1. Refresh the page.
2. Check your internet connection.
3. Try using a different browser (Chrome, Safari, Edge).
4. Clear your browser cache and cookies.
5. If the issue continues, try accessing the dashboard later.

8.1.2 Missing Data

Cause: The waist device may not have synced yet or the base station is out of range.

Solution:

1. Ensure the waist device is worn properly and turned on.
2. Make sure the base station is plugged in and within Bluetooth range.
3. Wait at least one hour to allow the next scheduled sync.
4. Refresh the dashboard afterward.

8.1.3 Inaccurate Scores

Cause: The device might not have been worn properly or didn't detect movement accurately.

Solution:

1. Double-check that the device is attached securely at the waist.
2. Avoid loose clothing that may interfere with motion detection.
3. Make sure the device was worn consistently throughout the day.
4. Contact your therapist if the data still appears off.

8.2 Device-Related Issues

8.2.1 Device Not Syncing

Cause: The base station may be unplugged, or Bluetooth range is too far.

Solution:

1. Make sure the base station is powered on.
2. Keep the waist device within 5–10 meters of the base station.

3. Try a manual sync (if your base station has a button for it).
4. Wait 15 minutes and check the dashboard again.

8.2.2 Battery Issues

Cause: The waist device battery may be low or dead.

Solution:

1. Recharge or replace the device battery as instructed.
2. Ensure the device is fully charged before daily use.
3. If the device does not power on after charging, contact support.

8.3 Who to Contact If the Issue Persists

If none of the steps above resolve the issue, please contact your therapist or the support team assigned to your care. They can:

- Remotely check device status
- Help with data corrections
- Arrange replacements or updates

9. Contact and Support

9.1 Helpdesk Contact

If you encounter any issues with the dashboard, device syncing, or your account, our helpdesk is here to assist you.

- Email: support@hipperacademy.nl

9.2 Support Hours

Our team is available on weekdays from 09:00 to 17:00 (CET).

Outside these hours, you're welcome to leave a message — we'll get back to you as soon as possible.

9.3 Response Time

We aim to respond to most support requests within 1–2 business days.

Critical issues, such as login problems or missing data, are prioritized and typically addressed within 24 hours.

9.4 Feedback Submission

Your feedback is always appreciated. If you'd like to report a bug, make a suggestion, or help us improve the platform:

- Send us a message at: info@hippertext.nl